

How to install the LORETA KEY software, version >2017


Prerequisites:

- 1. Windows 7 or higher (will not work correctly with WindowsXP or lower).**
- 2. The user must be logged into windows as an administrator with full rights as Administrator.**
- 3. The user must be familiar with the use of “Windows File Explorer”.**

Download:


- 1. Download the loreta key self-extracting file (“LoretaKey.exe”) to your computer.
Suggested locations: hard disk drive “C” or “D”, or to the “Desktop”.**
- 2. Open the program “Windows File Explorer”, and navigate to the location of “LoretaKey.exe”.**

Navigate here (in this example it is drive "D"):

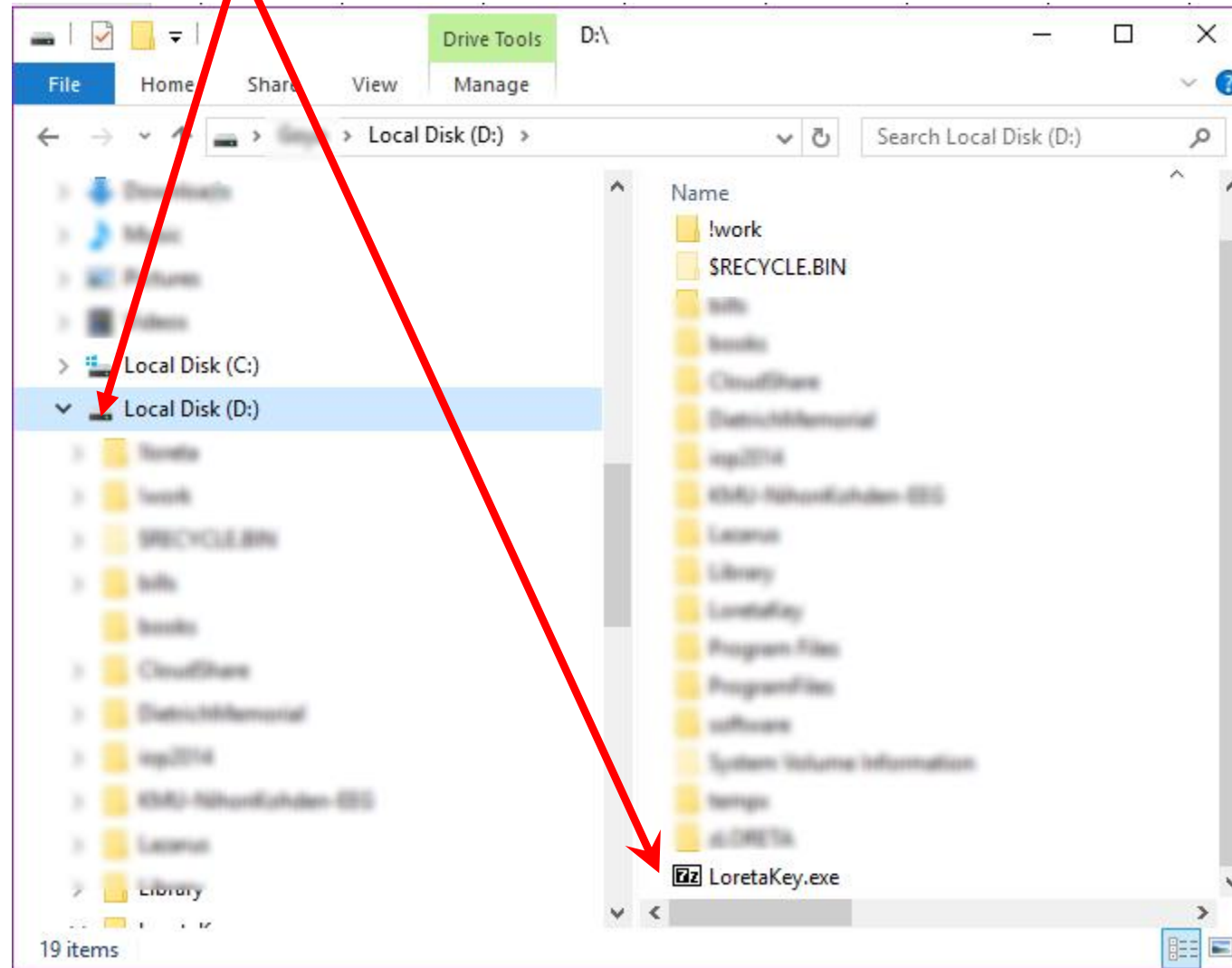


Right-click

Run as administrator

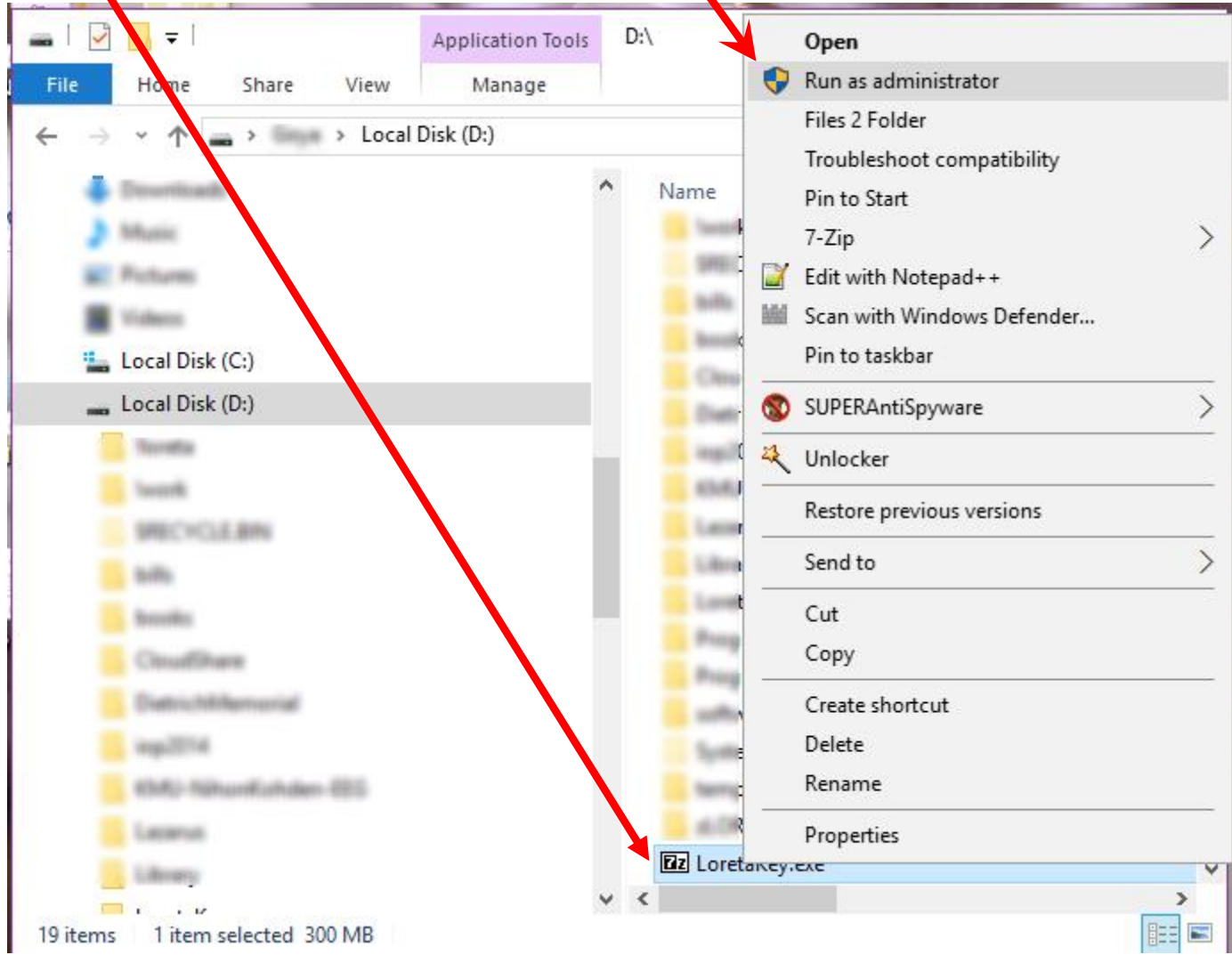

Click “EXTRACT”


Note: if you are an advanced user, you can decompress to a place of your own choosing. Be sure to not forget the location.

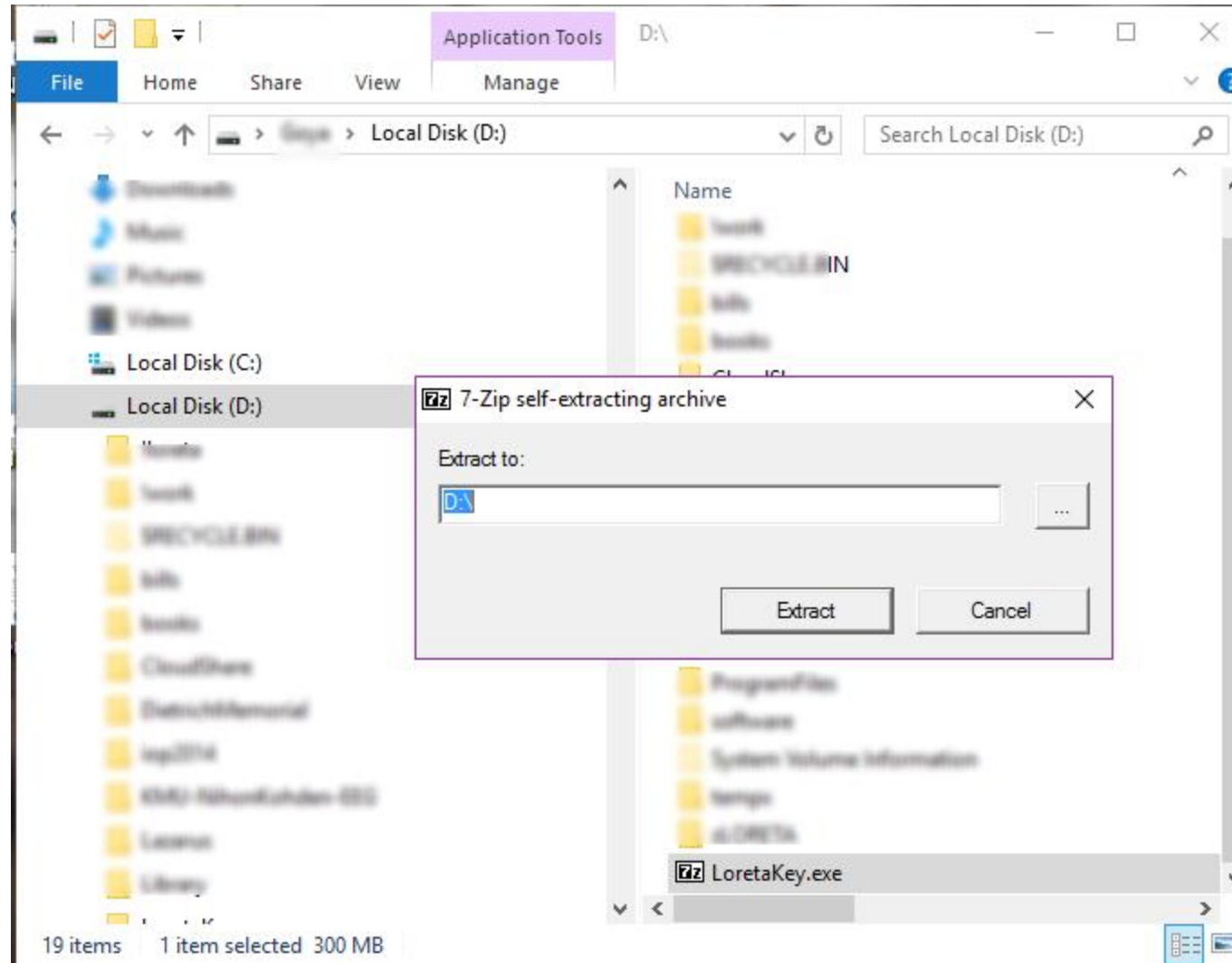


Select “Show password”


Type in the password (case sensitive): ZeroErrorLoreta

Click OK (the self extracting program starts to create the LoretaKey folder, please wait a bit...)

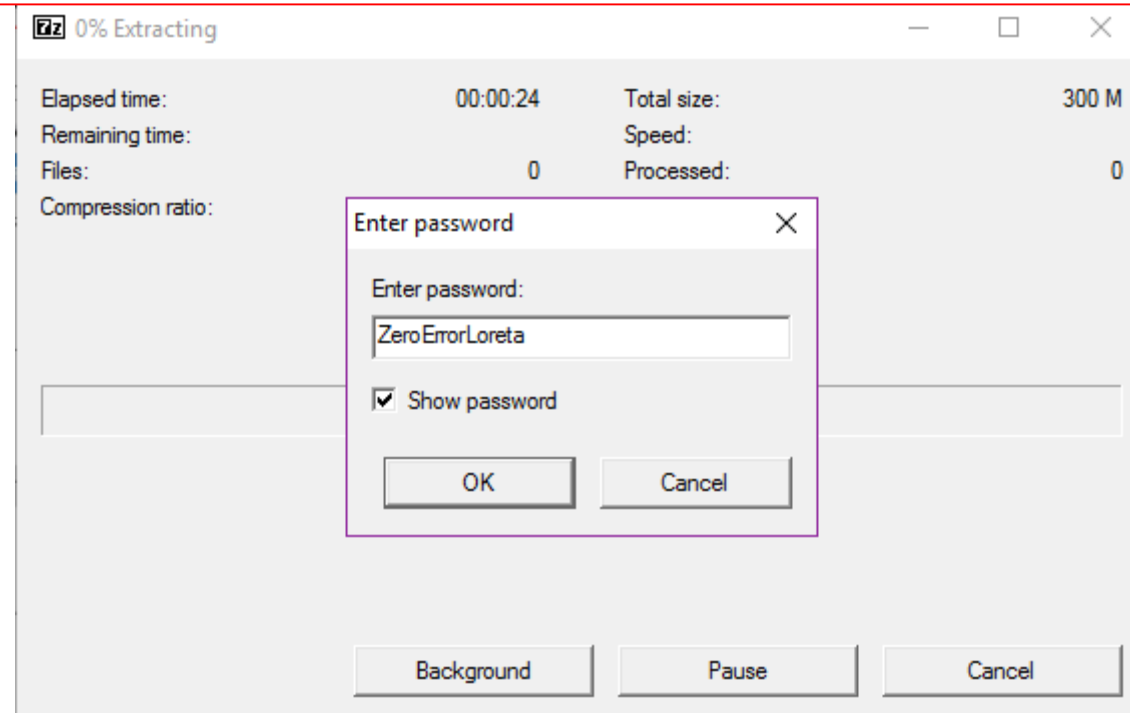


Note the newly created folder "LoretaKey" under drive "D" (in this example)

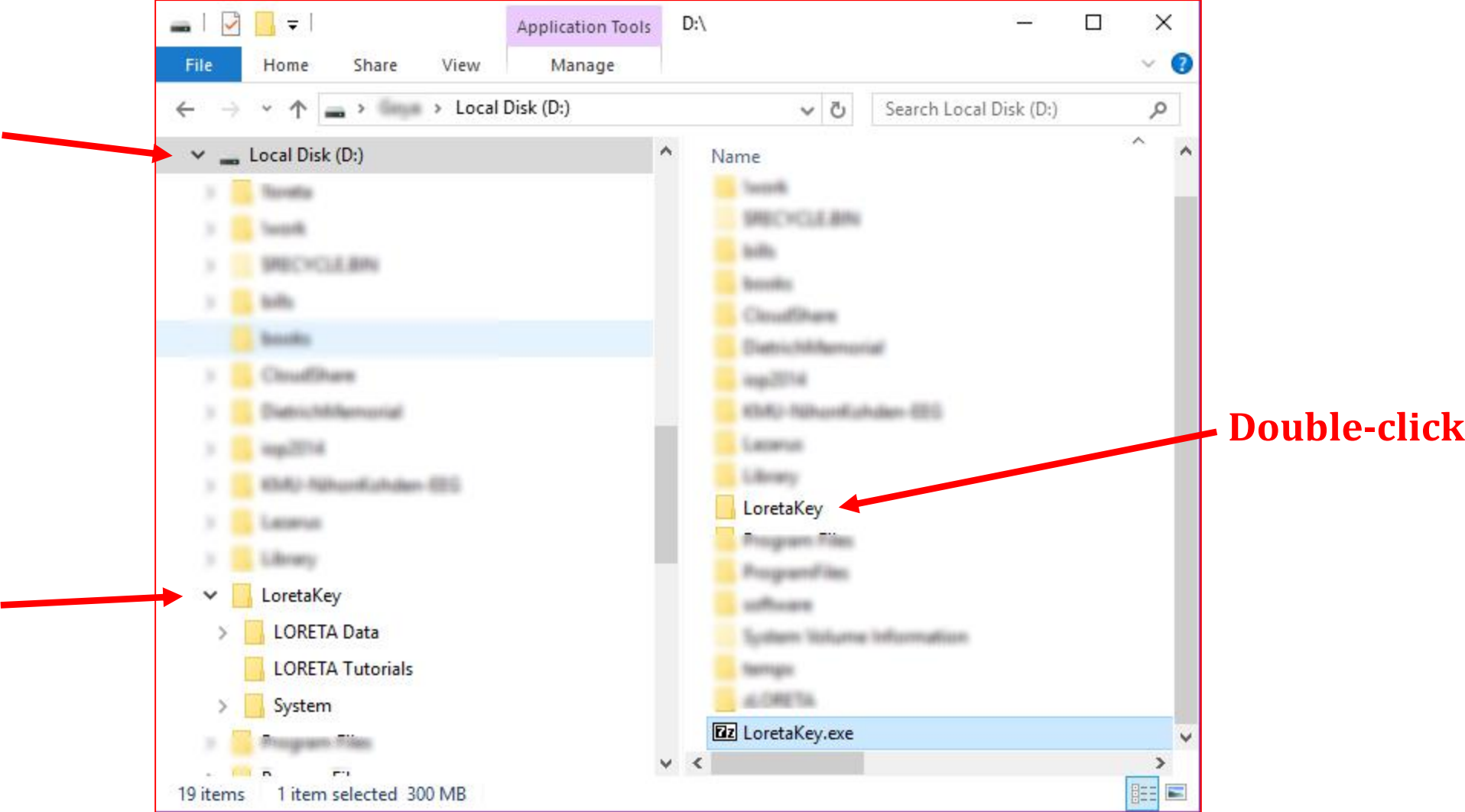

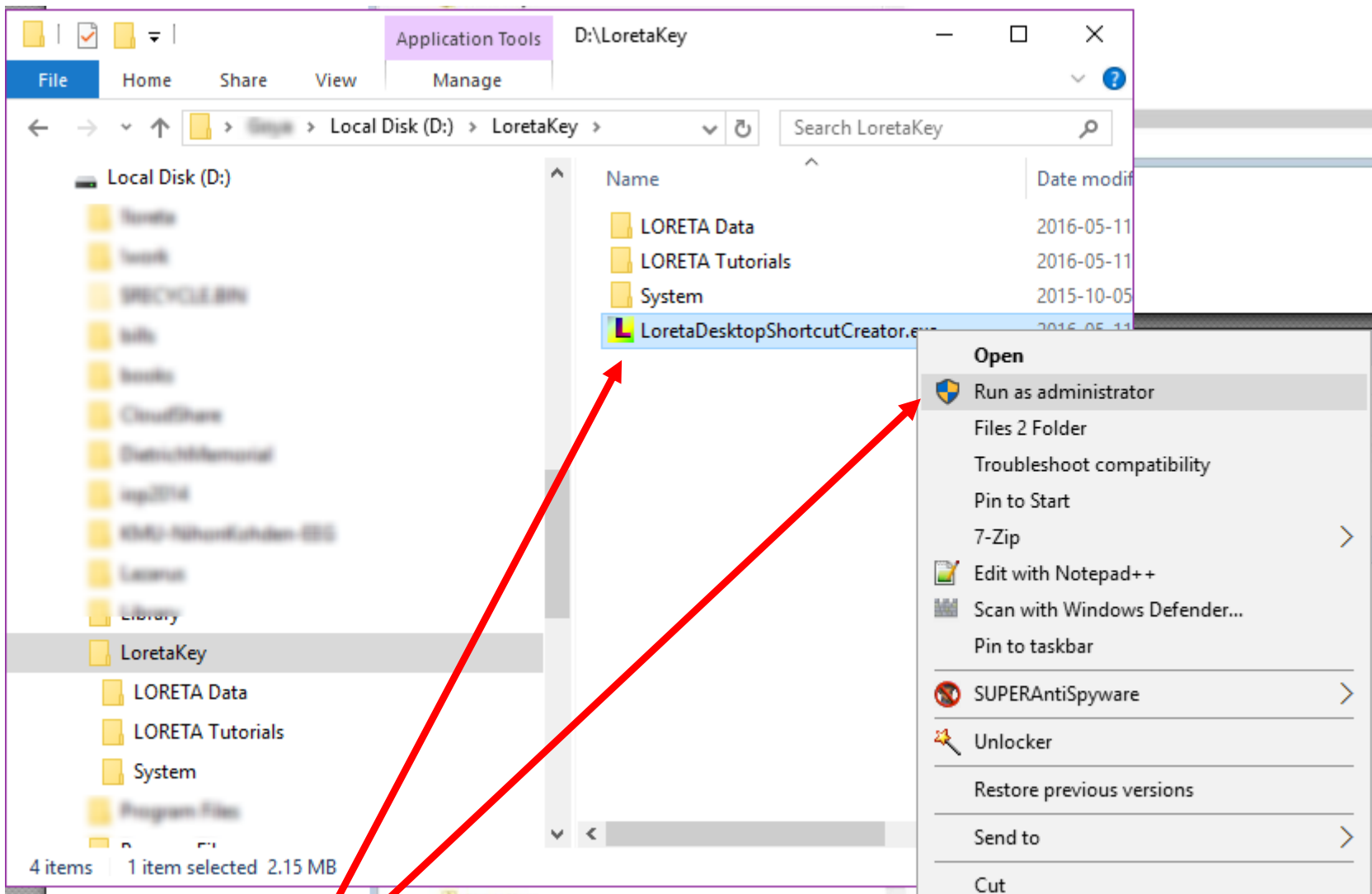


Right-click on the program “LoretaDesktopShortcutCreator.exe” and run as administrator

Click here

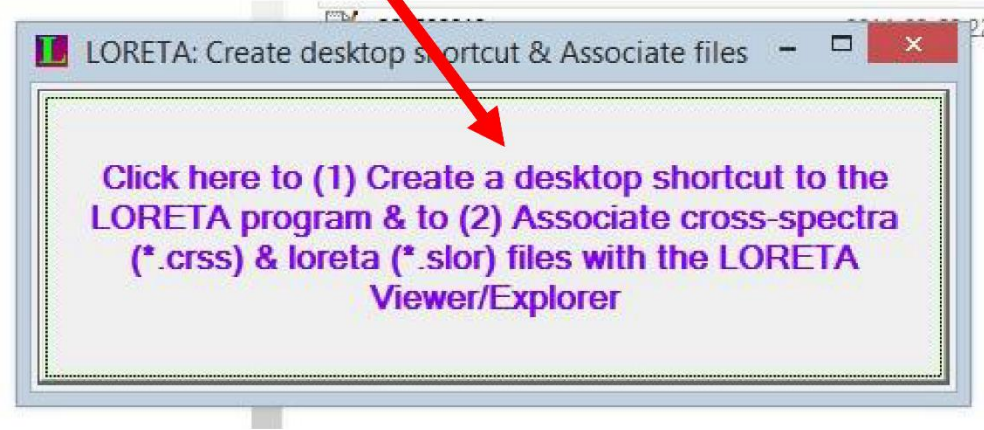

This will commence the process that first asks for the user's agreement with the loreta key licence, then creates a DESKTOP shortcut link to the loreta program, and then registers the file types "*.SLOR" and "*.CRSS" with the loreta viewer / explorer

This will commence the process that first asks for the user's agreement with the loreta key licence, then creates a DESKTOP shortcut link to the loreta program, and then registers the file types "*.SLOR" and "*.CRSS" with the loreta viewer / explorer

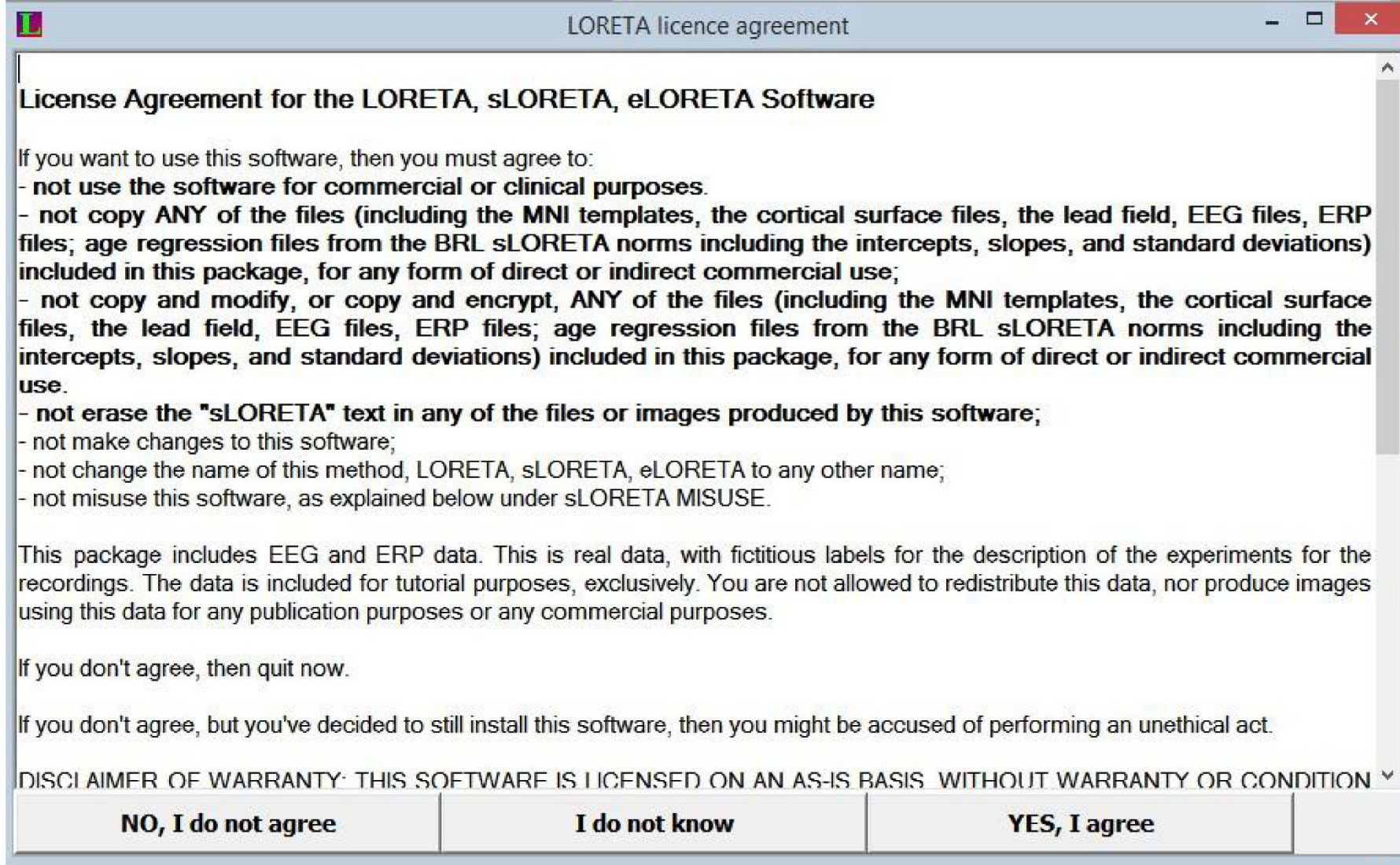

After agreement, click
ok to continue

3RL sLORETA norms including the intercept
m of direct or indirect commercial use;
d encrypt, ANY of the files (including the M
RP files; age regression files from the BI
viations)




ty of the this sc

ORETA, sL name;
elow unde

lata. This is real data, with fictitious labels for the
ial purposes, exclusively. You are not allowed to r
se or any commercial purposes


ie, for any form of direct or indirect commercial use;
r, or copy and encrypt, ANY of the files (including the MNI templates, the c
EEG files ERP files; age regression files from the BRL sLORETA norms
I stand

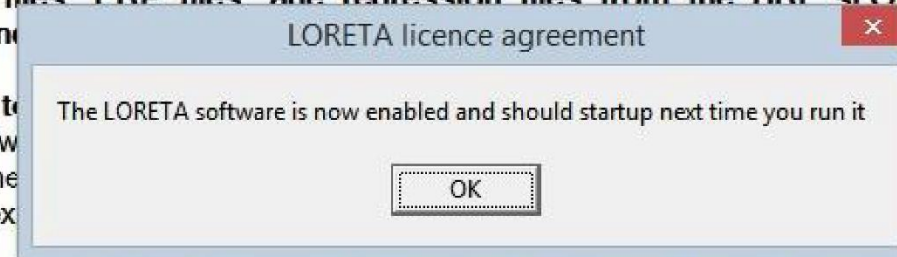


ETA* to
s softw
this me
, as ex

EG and ERP data. This is real data, with fictitious labels for the description of the exp
cluded for tutorial purposes, exclusively. You are not allowed to redistribute this data, nor


click ok to continue

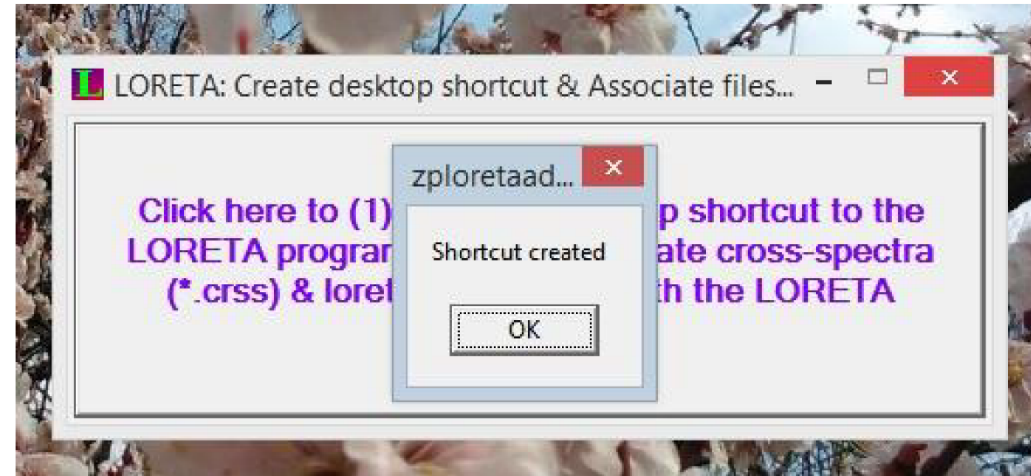


click ok to continue

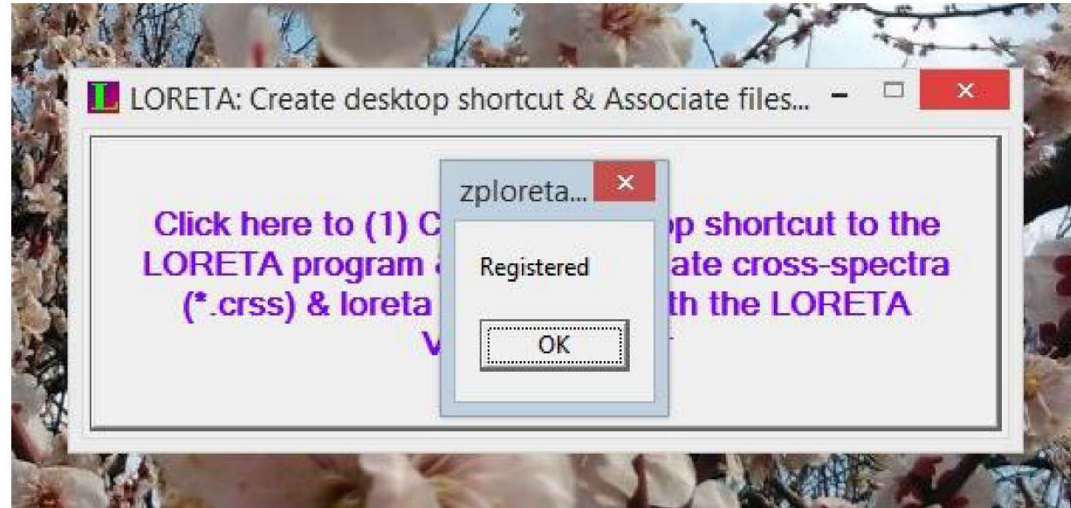

**Now your desktop will
have a shortcut link to
the loreta program**

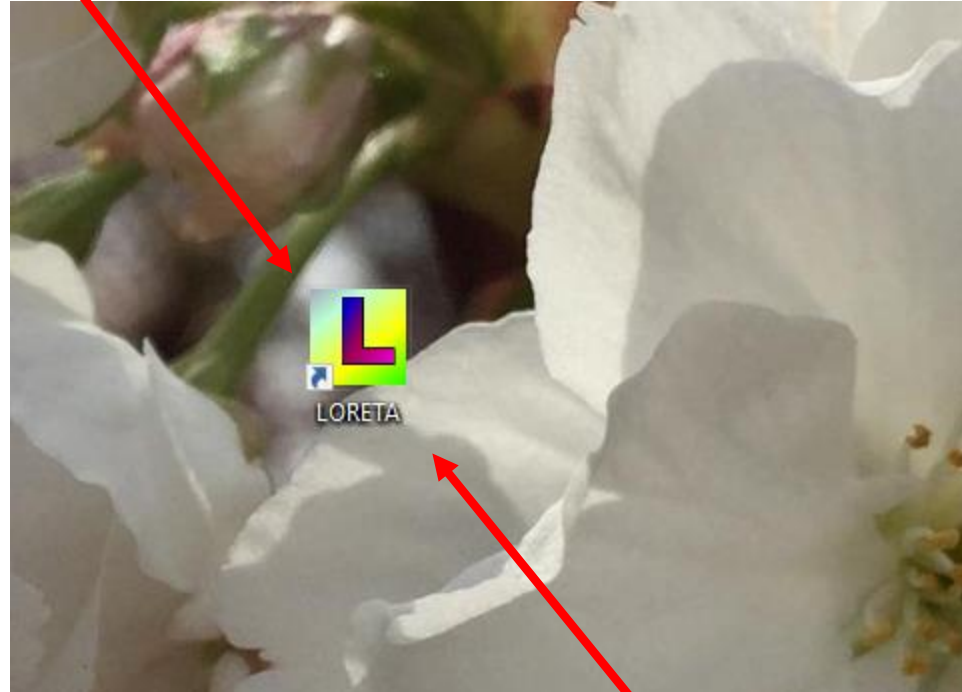

Double-click to run loreta

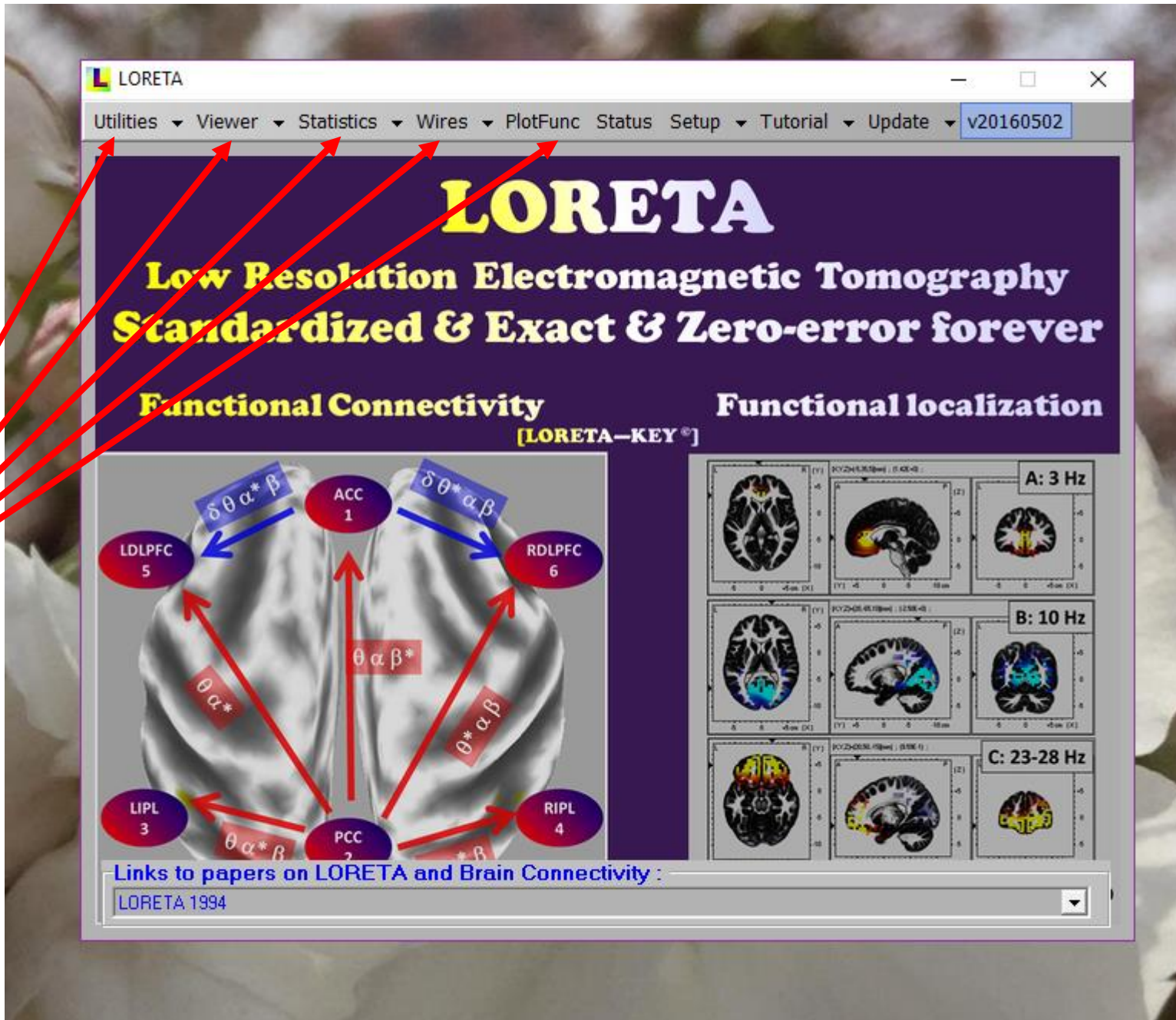

Now you are ready to use the program

Notes:

The folder “LORETA tutorials” contains the tutorials explaining how to use some selected parts of the many options of the software.

The folder “LORETA data” contains recordings and additional material, most of which is used in the tutorials.

